

Gold Anchor is ideal for breast biopsy marking with its unique and patented design. When deployed, Gold Anchor attaches immediately into the tissue.

Gold Anchor can be implanted percutaneously, guided by stereotactic mammography, ultrasound or CT.

Gold Anchor™



**breast
biopsy
marker**



Minimally invasive



Instant stability



Great visibility




Minimally Invasive

- Our thin needles easily fit needle guides used for biopsies.
- Ideal for placement in nodes.
- Easily done without the need for local anaesthesia.



20 G Gold Anchor

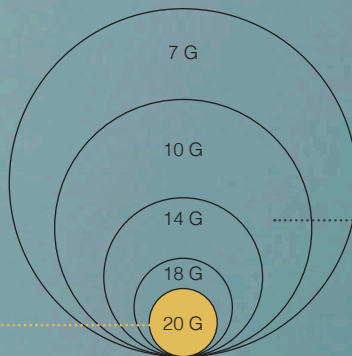


7 G-18 G Traditional breast biopsy markers

Instant Stability

Multiple cut-outs allow the marker to fold (Int'l patents). The marker is passive and will form different shapes depending on implantation technique.

Gold Anchor gets a great tissue attachment when the marker folds. The ball shaped marker becomes thicker than the needle tract.

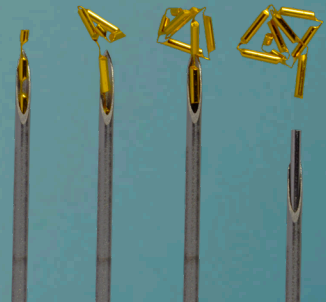

Great visibility

Gold Anchor is clearly visible on mammography (kV), CT, ultrasound, and MRI thanks to its unique marker material, a highly biocompatible alloy of pure gold and 0.5 or 1.5% pure iron (Int'l patents).

BM24.2



*Scan code
to request
samples*

Contact
info@goldanchormarker.com
order@goldanchormarker.com
www.GoldAnchorMarker.com

